


STERYFLUS TSE

- Easy integrity testable in situ
- Repeatedly steamable in situ or in autoclave
- Sanitizable
- Thermowelded construction
- EC-listed materials for Food contact
- FDA-listed materials per CFR21
- Bio-Safety per USP—Plastics

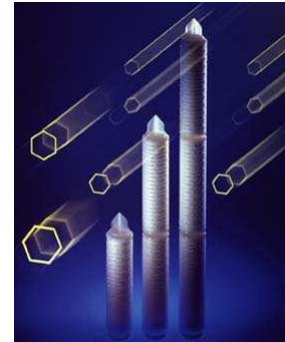

STERYFLUS TSE cartridge is final sterilizing filter for general application in pharmaceutical and food & beverage industries.

Steryflus cartridges have liquid absolute filtration ratings from 0,1 micron to 0,65 micron; 0,1 micron and 0,2 micron are available with single and double membrane layers.

Manufacturing is completed in a controlled environment; each cartridge is integrity tested and the limits of acceptability are monitored on regular basis by bacteria challenge test.

MATERIALS OF CONSTRUCTION

Filter media	Asymmetric PES membrane
Upstream supports	polyester
Downstream supports	polyester
Internal Core	polypropylene
External Cage	polypropylene
End caps / Adapters	polyester

FOOD-SAFETY

STERYFLUS TSE filter elements meet (EU) regulation 10/2011 and its subsequent amendments and regulations (EC) 1935/2004 and 1895/2005.

BIO-SAFETY

Filter media and components pass USP Biological Reactivity and Chemical-Physical tests for CLASS VI plastics.

QUALITY STANDARDS

Produced under a certified Quality System to guarantee traceability of manufacturing records and integrity testing results.



OPERATING CONDITIONS

- max. continuous temperature	65 °C
- max. cumulative time of steam sterilization	13 hours at 125 °C with cycles of 30 minutes
- sanitization with hot water	80 °C max
- sanitization with chemicals	can be sanitized by standard chemical agents
- max. differential pressure	5,0 bar at 25 °C - 1,7 bar at 80 °C
- recommended change out differential pressure	2,0 bar at 25 °C
- recommended rinse up volume	3 liters / 10" cartridge


CODE	ABSOLUTE FILTRATION RATING IN LIQUIDS	BACTERIAL RETENTION >10 ⁷ CFU/cm ² *	ACCEPTABLE LIMIT FOR DIFFUSION TEST WITH WATER FOR 10" CARTRIDGE (ml/min)
ZXX	0,1 µm double membrane	Acholeplasma laidlawii	≤ 20 @ 3,5 bar
ZX	0,1 µm	Acholeplasma laidlawii	≤ 20 @ 3,5 bar
ZYY	0,2 µm double membrane	Brevundimonas diminuta	≤ 16.5 @ 2,3 bar
ZY	0,2 µm	Brevundimonas diminuta	≤ 16.5 @ 2,3 bar
ZT	0,45 µm	Serratia marcescens	≤ 13 @ 1,7 bar
ZK	0,65 µm	N.A	≤ 20 @ 1,1 bar

* as per ASTM F838-05

WATER FLOW RATE FOR 10" CARTRIDGE

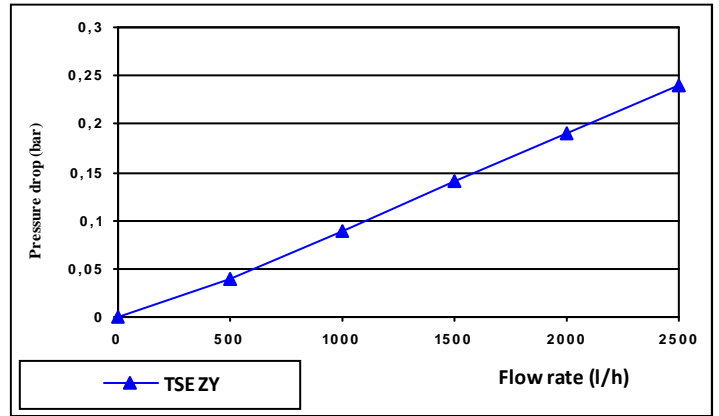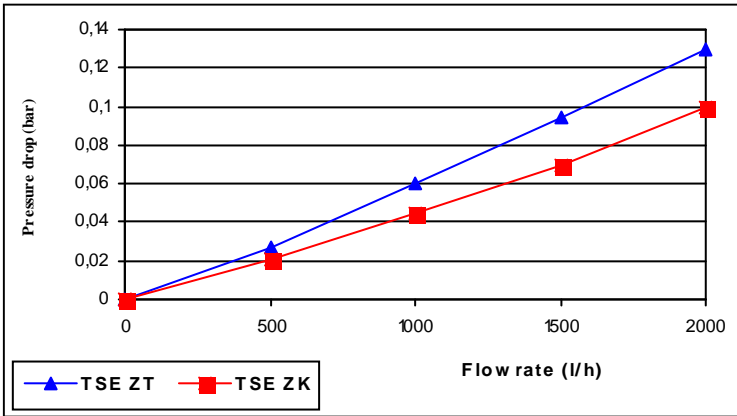

STERYFLUS TSE ORDERING INFORMATION

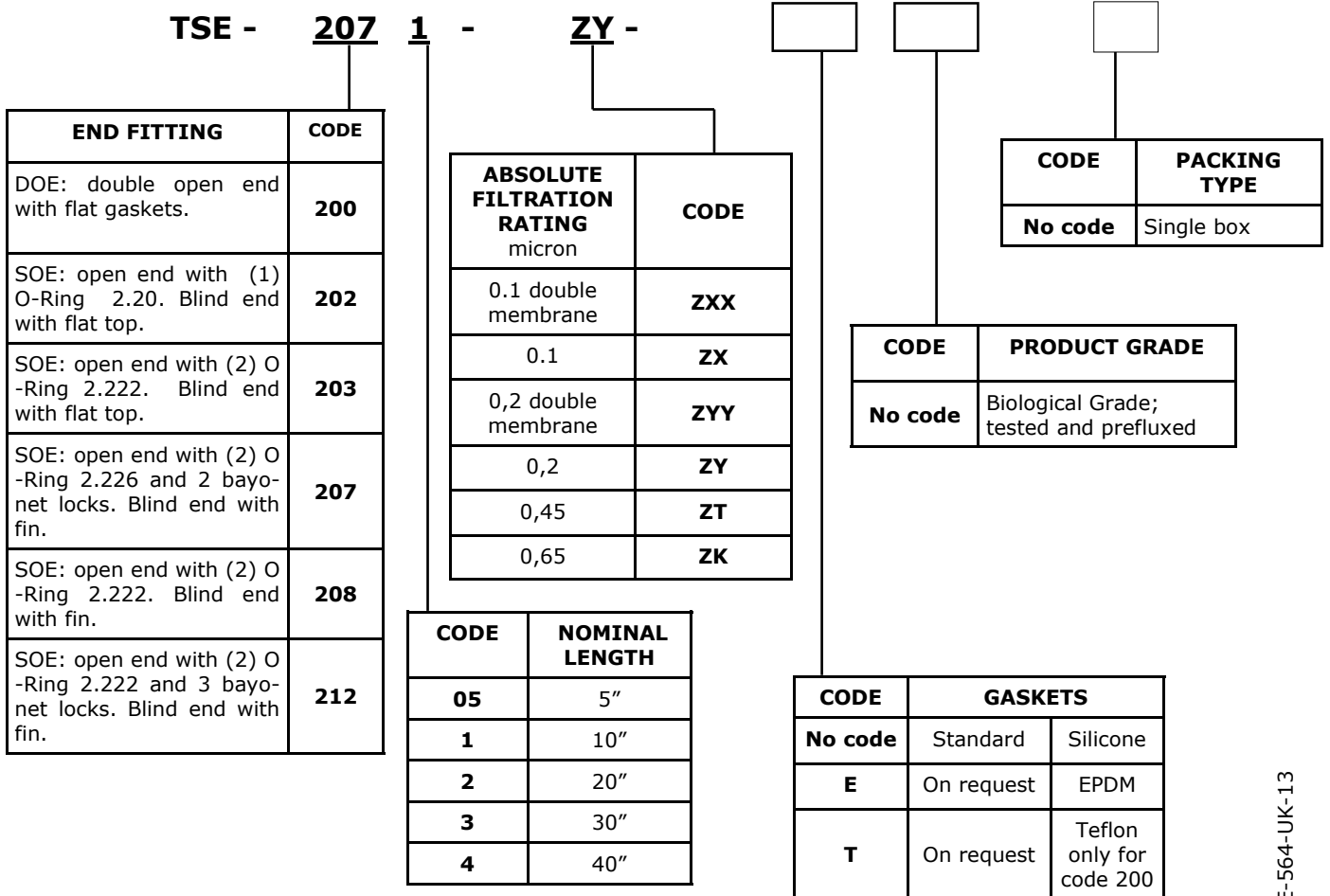

DS-TSE-564-UK-13

Data contained in this bulletin are informative and subject to change without notice.
 User is responsible for determining whether the product is fit for particular purpose and suitable for User's method of application.



Bea Technologies Spa Via Newton, 4 - 20016 Pero (Milano) ITALY
 Tel +39 02 339271 FAX +39 02 3390713 e-mail: info@bea-italy.com
 web: www.bea-italy.com